

If this is the first time logging in with this device, you must first complete the steps in the [First-Time Setup - MS Windows](#) guide.

Accessing GE PACS Web

Note: Access to GE PACS requires specific permissions. Please submit a Children's Service Desk **Account Request** for **Remote Access** to GE PACS, or call Children's IT Service Desk.

- 1 Log into Secure Access at <https://www.childrensmn.org/signin>
- 2 Click the link to **GE PACS Web**. 
- 3 Click the **link** for the **msxml.msi** install and then click **Run**.


Status	Check	Value	Description/Solution
ok	JavaScript	is enabled	
ok	Browser type/version	IE 7.0	
ok	Script Engine	11.0	
ok	Userdata persistence	is enabled	
✗	ActiveX	-	The installation of the following controls failed: - MSXML (please download the msxml.msi installer to install MSXML manually.) - Pictorial Index - Viewport

There is an error, Centricity Enterprise Web cannot go on.
Please refer to Chapter "Installing the Centricity Web Client" of the Installation Guide for further guidance or contact your system administrator.

[Click here to try again.](#)

Do you want to run or save **msxml.msi** (5.04 MB) from **gepacsweb.childrensmn.org**?


- 4 Continue through the installation and when prompted, click **Install Now**

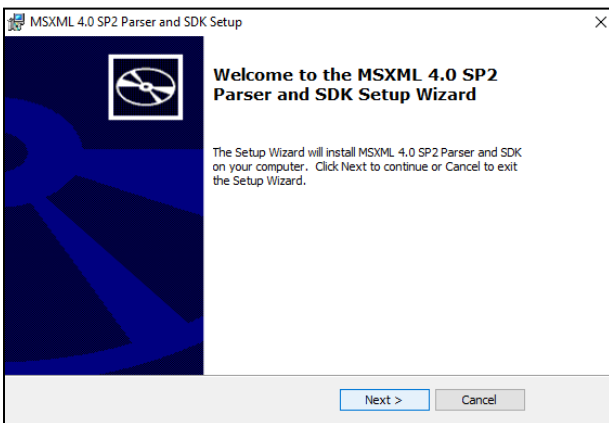

MSXML 4.0 SP2 Parser and SDK Setup

Welcome to the MSXML 4.0 SP2 Parser and SDK Setup Wizard

The Setup Wizard will install MSXML 4.0 SP2 Parser and SDK on your computer. Click Next to continue or Cancel to exit the Setup Wizard.

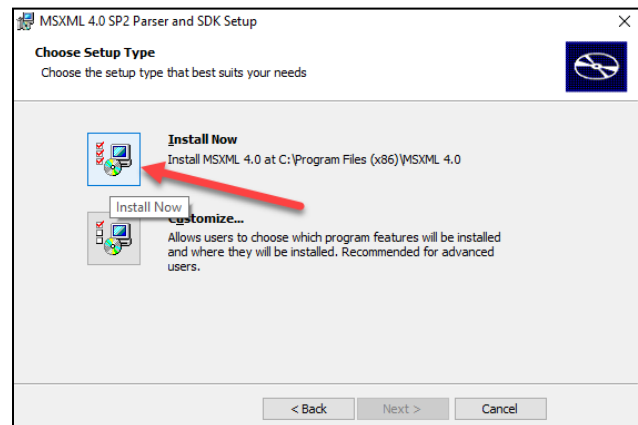

MSXML 4.0 SP2 Parser and SDK Setup

Choose Setup Type

Choose the setup type that best suits your needs

Install Now
Install MSXML 4.0 at C:\Program Files (x86)\MSXML 4.0


Install Now and Customize...
Allows users to choose which program features will be installed and where they will be installed. Recommended for advanced users.

- 5 Next, configure Internet Explorer to allow ActiveX, by going to **Internet Options** (gear symbol)




Internet Explorer window showing the Internet Options menu. The menu is open, and "Internet options" is highlighted. The background shows the GE Medical Systems Information Technologies Centricity Enterprise Web page.

6 Click the **Security tab**, highlight **Trusted Sites**, then click **Sites**.

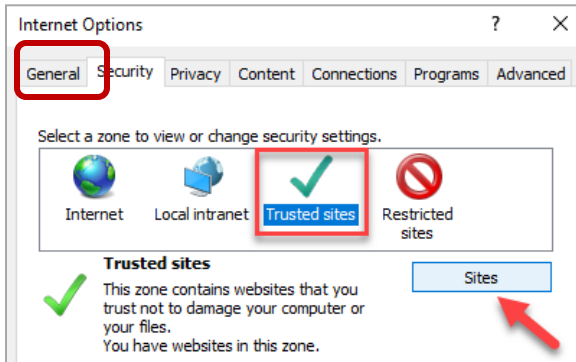

7 Uncheck “Require server verification...” and add <http://gepacsweb.childrensmn.org>

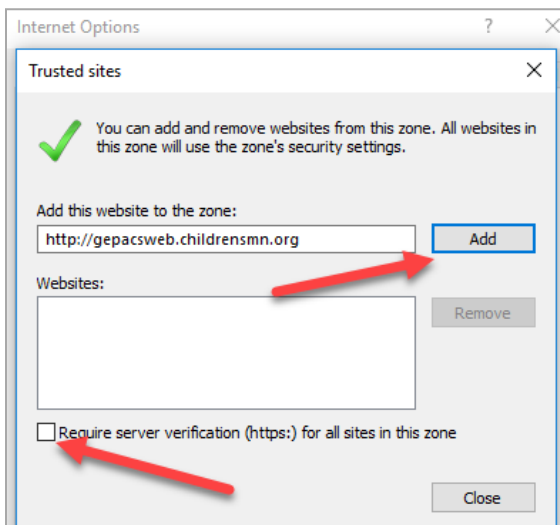

8 Now you can **install ActiveX**, and the **AMIVIEWER** and **PictorialIndex** add-on modules.

Status	Check	Value	Description/Solution
OK	JavaScript	is enabled	
OK	Browser type/version	IE 7.0	
OK	Script Engine	11.0	
OK	Userdata persistence	is enabled	
✗	ActiveX	-	The installation of the following controls failed: - Pictorial Index - Viewport

There is an error, Centricity Enterprise Web cannot go on.

Please refer to Chapter "Installing the Centricity Web Client" of the Installation Guide for further guidance or contact your system administrator.

[Click here to try again.](#)

Internet Explorer blocked this website from installing an ActiveX control. [What's the risk?](#) Install

This webpage wants to run the following add-on: 'AMIVIEWER ActiveX Control Module' from 'GE Medical Systems (unverified publisher)'. [What's the risk?](#) Allow

This webpage wants to run the following add-on: 'PictorialIndex ActiveX Control Module' from 'GE Medical Systems (unverified publisher)'. [What's the risk?](#) Allow

Click the Centricity **Try Again** link and sign into GE PACS with your Children's MN username /password.